

# Reliable Collection and Transport of Samples Ensures Accurate Patient Test Results.



## Media

The media tube is designed with a spiral groove, round bottom media chamber which keeps the sample immersed and protected in the media while allowing easy racking. The color coded caps are easily identified and feature a protective biohazard overskirt. Each tube has a patient identification label and product code, lot number and expiry date.

A full range of media formulations include: Modified Liquid Stuart's, Modified Amies Clear, Modified Amies Charcoal, Modified Liquid Amies.

## Modified Liquid Stuart's and Modified Liquid Amies

Provides excellent recovery in routine bacteriology and is compatible with gram stain and many rapid screening methods. Shelf life is 18 months.

## Modified Amies Clear Gel

Provides excellent survival for a wide range of bacteria and is compatible with gram stain. Shelf life is 24 months.

## Modified Amies Charcoal Gel

Further enhances the survival of fastidious organisms including Neisseria gonorrhoeae by neutralizing any fatty acids that may be present in the specimen. Shelf life is 24 months.



[www.starplexscientific.com](http://www.starplexscientific.com)  
1-800-665-0954

# Reliable Collection and Transport of Samples Ensures Accurate Patient Test Results.

## Starswabs II®

The media can be paired with a variety of swab tips and shafts. Studies available upon request.

- Single swab, plastic shaft with rayon tip
- Double swab, plastic shaft with rayon tip
- Straight Wire shaft with rayon tip

## Packaging

Sterile peel pouch single or double packed swabs are sealed into a metalized bag.

- Premium thickness bag creates a barrier environment to prevent premature medium dehydration
- Bag also features a Ziplock closure which promotes resealing the bag to maintain swabs at peak performance
- Metalized bag protects media from harmful UV rays


## Dry Swabs

Also available are single and double rayon tipped swabs, without media, in the spiral, round bottom tube with a white plug cap. Sterile with patient identification label. Used in forensic collections.

## Anaerobic Transport System

*Especially for fastidious anaerobes*

The System has an indicator that will turn lavender if the anaerobic condition is compromised. Furthermore, letting oxygen into the tube can be avoided by injecting fluids through the rubber septum of the screw cap with a needle and syringe. Both a 7ml and 5ml media fill are available. Two sterile scored swabs packaged with each tube. 18 month shelf life.

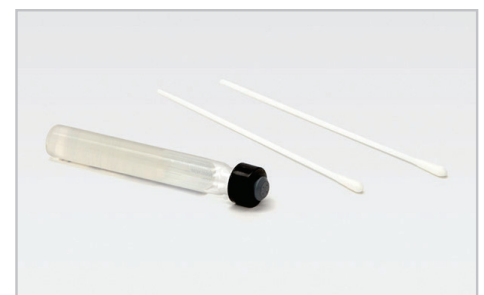

[www.starplexscientific.com](http://www.starplexscientific.com)  
1-800-665-0954